







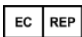


NATRACURE™ ADVANCED COOLING GEL MITTS 1 Pair 11.8" x 7.9" x 0.8"	REF 761-00NC1
NATRACURE™ AVANCEREDE KØLENDE GELVANTER 1 par 29,97 x 20,07 x 2,03 cm	GANTS DE GEL DE REFROIDISSEMENT NATRACURE™ ADVANCES 1 Paire 29,97 x 20,07 x 2,03 cm
NATRACURE™ HANDSCHUHE IN FORTSCHRITTLICHER KÜHLGELTECHNOLOGIE 1 Paar 29,97 x 20,07 x 2,03 cm	NATRACURE™ QUANTI IN GEL A RAFFREDDAMENTO AVANZATO 1 Paio 29,97 x 20,07 x 2,03 cm
NATRACURE™ AVANSERTE KJØLEGEHANSKER 1 par 29,97 x 20,07 x 2,03 cm	NATRACURE™ LUVAS DE GEL DE REFRIGERAÇÃO AVANÇADA 1 Par 29,97 x 20,07 x 2,03 cm
MITONES DE GEL FRIO AVANZADAS NATRACURE™ 1 par 29,97 x 20,07 x 2,03 cm	NATRACURE™ VANTAR MED AVANCERAD KYLANDE GEL 1 par 29,97 x 20,07 x 2,03 cm

 QualRep Services BV Utrechtseweg 310 – B42, 6812 AR, Arnhem – the Netherlands Phone: +31 20 788 2630 globalreg@qservegroup.com	 PolyGel, LLC 1 Danforth Drive, Easton, PA 18045 USA	
QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Holland Telefon: +31 20788 2630	QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Pays-Bas Téléphone: +31 20788 2630	 YYYY-MM-DD
QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Paesi Bassi Telefono: +31 20788 2630	QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Niederlande Telefon: +31 20 788 2630	    
QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Holanda Telefone: +31 20 788 2630	QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Nederland Telefon: +31 20788 2630	UDI BAR CODE (FPO)
QualRep Services BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Nederlanderna Telefon: +31 20788 2630	Servicios de QualRep BV Utrechtseweg 310 - B42, 6812 AR, Arnhem - Países Bajos Teléfono: +31 20788 2630	

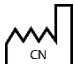 Officiële geautoriseerde distributeur E.B.Materials b.v. G. Gezellelaan 36/2 3550 Heusden-Zolder België	www.belgomedical.com/store-1/ Tel +32 11 53 77 47
--	---

NATRACURE™ ADVANCED COOLING GEL MITTS**REF 761-00NC1**

1 Pair | 11.8" x 7.9" x 0.8"

Designed to help:

- Soothe hot, painful fingers during and after chemotherapy
- Reduce neuropathy pain and nail loss*

Key features:

- New GEL/ICE™ Formula for long cold therapy treatments
- Wide opening allows for easy on-and-off
- Convenient hook-and-loop enclosure at wrist
- Durable, stain-resistant exterior
- Soft, inner fabric

Intended Use(s): Provide cold therapy to specific body parts to support healing post-injury, surgery and/or chronic pain.

Intended User(s): Intended for use by individuals with neuropathy pain and discomfort in hands and fingers.

Directions: Place both mitts in the refrigerator or freezer, either wrapped in a paper towel or placed inside the accompanying plastic cover. Chill for a minimum of 2 hours. Place hands inside mitts comfortably and use hook-and-loop closure to secure the fit. Keep on for 20 minutes at a time.

Care: Outer and inner covers can be hand washed with soap and water. Air dry or pat dry until completely dry. Washing with soaps that have less scent and other possible additives minimize any possible reactions. Dispose of in accordance with local regulations.

Contraindication(s): Individuals using pacemakers or insulin pumps should NEVER use these products. Do not use during pregnancy. Do not apply cold to area of lost sensation. Do not use where there is poor circulation.

Warning(s):

1. FOR COLD USE ONLY!
2. Do not place product on an open wound.
3. If product causes skin irritation, discomfort or poor circulation occurs, discontinue use and notify your physician.
4. If you have diabetes, neuropathy, or poor circulation, consult a medical professional before using.
5. DO NOT STEP ON PACKS! Excessive weight may cause pack(s) to leak or rupture.
6. USE OVER CLOTHING OR PLACE FABRIC INTERFACE BETWEEN THE COLD PACK AND THE SKIN, PRIOR TO USING COLD THERAPY. POSSIBLE INTERFACES INCLUDE: A light towel, cloth, sheet/pillowcase (cut to size), 6 layers of paper towel, or Stockinette. Either wrap the interface around the hot/cold pack prior to attaching to the support (if included) or hold the interface on the skin while placing the product on the body.

Caution(s):

1. Consult your physician before using this product.
2. Only use this product for 20 minutes at a time on the affected area.
3. If at any time soreness, rash, or an increase in pain occurs, discontinue use and consult your physician.
4. Only use this product for its intended purposes as described in the instructions for use.
5. Do not use this product while using transdermal drug delivery skin patches, skin creams, balms, liniments or lotions.
6. If removable pack has been torn or has been compromised (and is leaking) discard immediately.
7. Always use a fabric interface.
8. Keep out of reach of children and pets.

Materials: PE/Nylon Loop Velvet, PTFE-coated nylon, Proprietary Gel Formula

UMDNS: 25806

Made in China
©2024 PolyGel®, LLC
PolyGel, LLC
1 Danforth Drive, Easton, PA 18045 USA
NATRACURE™ is a brand of PolyGel, LLC

NatraCure.com

* Studies that cite cold therapy for neuropathy pain and nail loss
<https://pubmed.ncbi.nlm.nih.gov/19171552/>
<https://pubmed.ncbi.nlm.nih.gov/32036471/>
<https://www.breastcancer.org/research-news/cold-therapy-helps-hair-skin-nail-effects>
<https://jamanetwork.com/journals/jamadermatology/article-abstract/2711245>

EC REP

QualRep Services BV
Utrechtseweg 310 – B42, 6812 AR, Arnhem – the Netherlands
Phone: +31 20 788 2630
globalreg@qservegroup.com

Officiële geautoriseerde distributeur
E.B.Materials b.v.
G. Gezellelaan 36/2
3550 Heusden-Zolder
België

www.belgomedical.com/store-1/ DANISH | DANSK
Tel +32 11 53 77 47

Officiële geautoriseerde distributeur
E.B.Materials b.v.
G. Gezellelaan 36/2
3550 Heusden-Zolder
België

www.belgomedical.com/store-1/
Tel +32 11 53 77 47

Not Made with Natural Rubber Latex



